


Use large block display to compare Bur identification system when using the schematic below for proper Bur usage

● VT5 Set    ○ VT8 Set    ● VS8 Set

Densifying Mode CCW (800-1500) RPMs / Cutting Mode CW (800-1500) RPMs																	
Astra Tech™			OsseoSpeed™ TX														
			Soft Bone							Hard Bone (Mandible)							
										In densifying mode make sure your osteotomy is 1.0 mm deeper than the actual implant final length. In extreme hard bone, utilize DAC (Densify After Cut) Protocol. <a href="#">Click here</a> to view PDF.							
Geometry	Major Ø	Minor Ø	Pilot	Bur 1	Bur 2	Bur 3	Bur 4	Bur 5	Densah® Bur Block Display	Pilot	Bur 1	Bur 2	Bur 3	Bur 4	Bur 5	Bur 6	Densah® Bur Block Display
Straight	3.0	3.0	Pilot	VT1828 (2.3)	VS2228* (2.5)	—	—	—		Pilot	VT1828 (2.3)	VT2535* (3.0)	—	—	—	—	
Straight	3.5	3.5	Pilot	VT1525 (2.0)	VT2535* (3.0)	—	—	—		Pilot	VT1525 (2.0)	VT1828 (2.3)	VT2535 (3.0)	VT2838* (3.3)	—	—	
Straight	4.0	4.0	Pilot	VT1828 (2.3)	VT2838 (3.3)	VS3238* (3.5)	—	—		Pilot	VT1828 (2.3)	VT2535 (3.0)	VT2838 (3.3)	VT3545* (4.0)	—	—	
Conical	4.5	3.5	Pilot	VT1525 (2.0)	VT2535 (3.0)	VT2838* (3.3)	—	—		Pilot	VT1525 (2.0)	VT2535 (3.0)	VT2838 (3.3)	VT3545* (4.0)	—	—	
Conical	5.0	4.0	Pilot	VT1828 (2.3)	VT2838 (3.3)	VT3545* (4.0)	—	—		Pilot	VT1828 (2.3)	VT2535 (3.0)	VT2838 (3.3)	VT3545 (4.0)	VT3848* (4.3)	—	
Straight	5.0	5.0	Pilot	VT1828 (2.3)	VT2838 (3.3)	VT3848 (4.3)	VS4248* (4.5)	—		Pilot	VT1828 (2.3)	VT2838 (3.3)	VT3545 (4.0)	VT3848 (4.3)	VT4555* (5.0)	—	

\*Denotes implant placement.

NOTE: Surgeon preference overrides this suggestive protocol

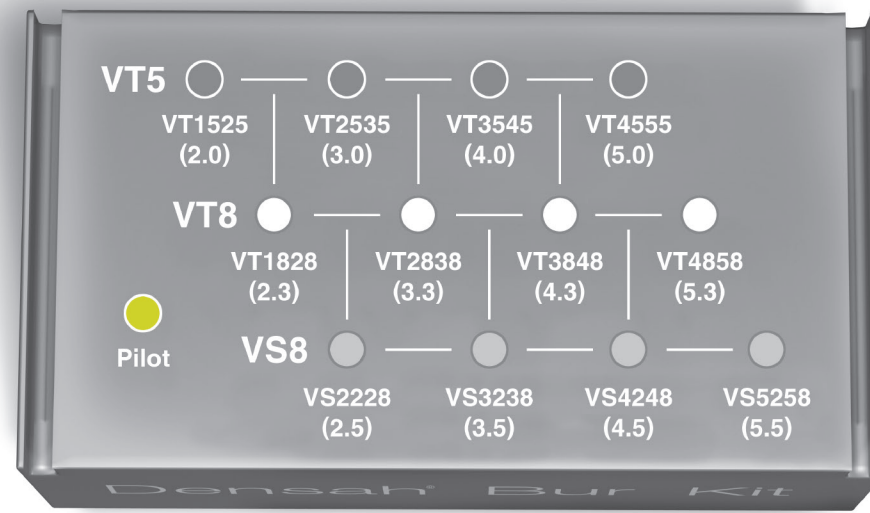

Use large block display to compare Bur identification system when using the schematic below for proper Bur usage

● VT5 Set    ○ VT8 Set    ● VS8 Set

Densifying Mode CCW (800-1500) RPMs / Cutting Mode CW (800-1500) RPMs																	
Astra Tech™			OsseoSpeed™ EV														
			Soft Bone							Hard Bone (Mandible)							
										In densifying mode make sure your osteotomy is 1.0 mm deeper than the actual implant final length. In extreme hard bone, utilize DAC (Densify After Cut) Protocol. <a href="#">Click here</a> to view PDF.							
Geometry	Major Ø	Minor Ø	Pilot	Bur 1	Bur 2	Bur 3	Bur 4	Bur 5	Densah® Bur Block Display	Pilot	Bur 1	Bur 2	Bur 3	Bur 4	Bur 5	Bur 6	Densah® Bur Block Display
Straight	3.0	3.0	Pilot	VT1828 (2.3)	VS2228* (2.5)	—	—	—		Pilot	VT1525 (2.0)	VT1828 (2.3)	VT2535* (3.0)	—	—	—	
Straight	3.6	3.6	Pilot	VT1525 (2.0)	VT2535 (3.0)	VT2838* (3.3)	—	—		Pilot	VT1525 (2.0)	VT1828 (2.3)	VT2535 (3.0)	VS3238* (3.5)	—	—	
Straight	4.2	4.2	Pilot	VT1828 (2.3)	VT2838 (3.3)	VS3238* (3.5)	—	—		Pilot	VT1828 (2.3)	VT2535 (3.0)	VT2838 (3.3)	VT3545* (4.0)	—	—	
Conical	4.2	3.6	Pilot	VT1828 (2.3)	VT2838 (3.3)	VS3238* (3.5)	—	—		Pilot	VT1828 (2.3)	VT2535 (3.0)	VT2838 (3.3)	VS3238* (3.5)	—	—	
Conical	4.8	4.2	Pilot	VT1525 (2.0)	VT2535 (3.0)	VT3545* (4.0)	—	—		Pilot	VT1525 (2.0)	VT2535 (3.0)	VT2838 (3.3)	VT3545* (4.0)	VT3848* (4.3)	—	
Straight	4.8	4.8	Pilot	VT1525 (2.0)	VT2535 (3.0)	VT3545 (4.0)	VT3848* (4.3)	—		Pilot	VT1525 (2.0)	VT2535 (3.0)	VT2838 (3.3)	VT3545 (4.0)	VT3848 (4.3)	VS4248* (4.5)	
Straight	5.4	5.4	Pilot	VT1828 (2.3)	VT2838 (3.3)	VT3848 (4.3)	VT4555* (5.0)	—		Pilot	VT1828 (2.3)	VT2838 (3.3)	VT3545 (4.0)	VT3848 (4.3)	VT4555 (5.0)	VT4858* (5.3)	

\*Denotes implant placement.

NOTE: Surgeon preference overrides this suggestive protocol